

Watmind SARS-CoV-2 antigen rapid test

Safe and practical
point-of-care assay



Reliable and cost-effective SARS-CoV-2 antigen rapid test

Ideal solution for fast and qualitative results in patients with onset of symptoms between 0 – 7 days.


Key benefits

- ✓ For professional use at point of care
- ✓ All sample collection and testing material included
- ✓ Saliva, naso-/oropharyngeal sampling
- ✓ Quick result after 10 minutes
- ✓ Single one-step buffer and mixing tubes
- ✓ Optional QC material available
- ✓ Based on colloidal gold immunoassay technology



Specifications*

Test	
Technology	Lateral flow sandwich immunoassay
Sample type	Saliva, naso-/oropharyngeal swab
Buffer volume	26 x 500 µL (ready-to-use buffer and mixing tube)
Sample volume	100 µL (3 drops)/test
Sensitivity	Naso-/oropharyngeal sample 87.0% (RT-PCR CT value ≤ 36); 95.15% (≤ 7 days symptom onset) Saliva sample 90.45% (RT-PCR CT value ≤ 36); 95.51% (≤ 7 days symptom onset)
Specificity	Naso-/oropharyngeal sample: 99.12% Saliva sample: 99.45%
Packaging specification	25 tests/kit

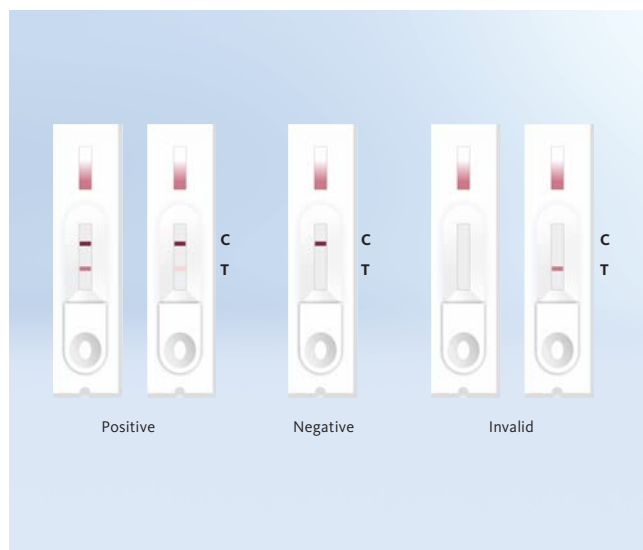

* For more details please refer to the IFU.

Distributor EMEA: Sysmex Europe GmbH

Bornbarch 1, 22848 Norderstedt, Germany · Phone +49 40 52726-0 · Fax +49 40 52726-100 · info@sysmex-europe.com · www.sysmex-europe.com

Authorised representative (EC REP): Shanghai International Holding Corp. GmbH (Europe)

Eiffestrasse 80, 20537, Hamburg, Germany

Manufacturer: Shenzhen Watmind Medical Co., Ltd.

Room 108, Shenzhen IC Design & Application Industrial Park Nr 1089, Chaguang Road, Nanshan District, Shenzhen, China

You will find your local Sysmex representative's address under www.sysmex-europe.com/contacts